



Swiss Institute of
Bioinformatics

SINGLE-CELL TRANSCRIPTOMICS WITH R

Single-cell RNA-seq and Cell Ranger

Deepak Tanwar
June 16-18, 2026

Why single-cell RNA-seq?

Why single-cell RNA-seq?

Single-cell RNA-seq (scRNA-seq) allows us to evaluate the transcriptome at the level of individual cells.



Why scRNAseq?

- To explore which cell types are present in a tissue

Why scRNAseq?

- To explore which cell types are present in a tissue
- To identify unknown/rare cell types or states

Why scRNAseq?

- To explore which cell types are present in a tissue
- To identify unknown/rare cell types or states
- To elucidate the changes in gene expression during differentiation processes or across time or states


Why scRNAseq?

- To explore which cell types are present in a tissue
- To identify unknown/rare cell types or states
- To elucidate the changes in gene expression during differentiation processes or across time or states
- To identify genes that are
 - differentially expressed in particular cell types
 - between conditions (e.g. treatment or disease)

Bulk vs single-cell

Bulk RNA-seq provides an overview of **average differences in gene expression**. For certain scenarios this has proven to be sufficient (i.e. biomarkers for cancer).

Bulk RNA-seq is useful if you are not expecting or not concerned about cellular heterogeneity

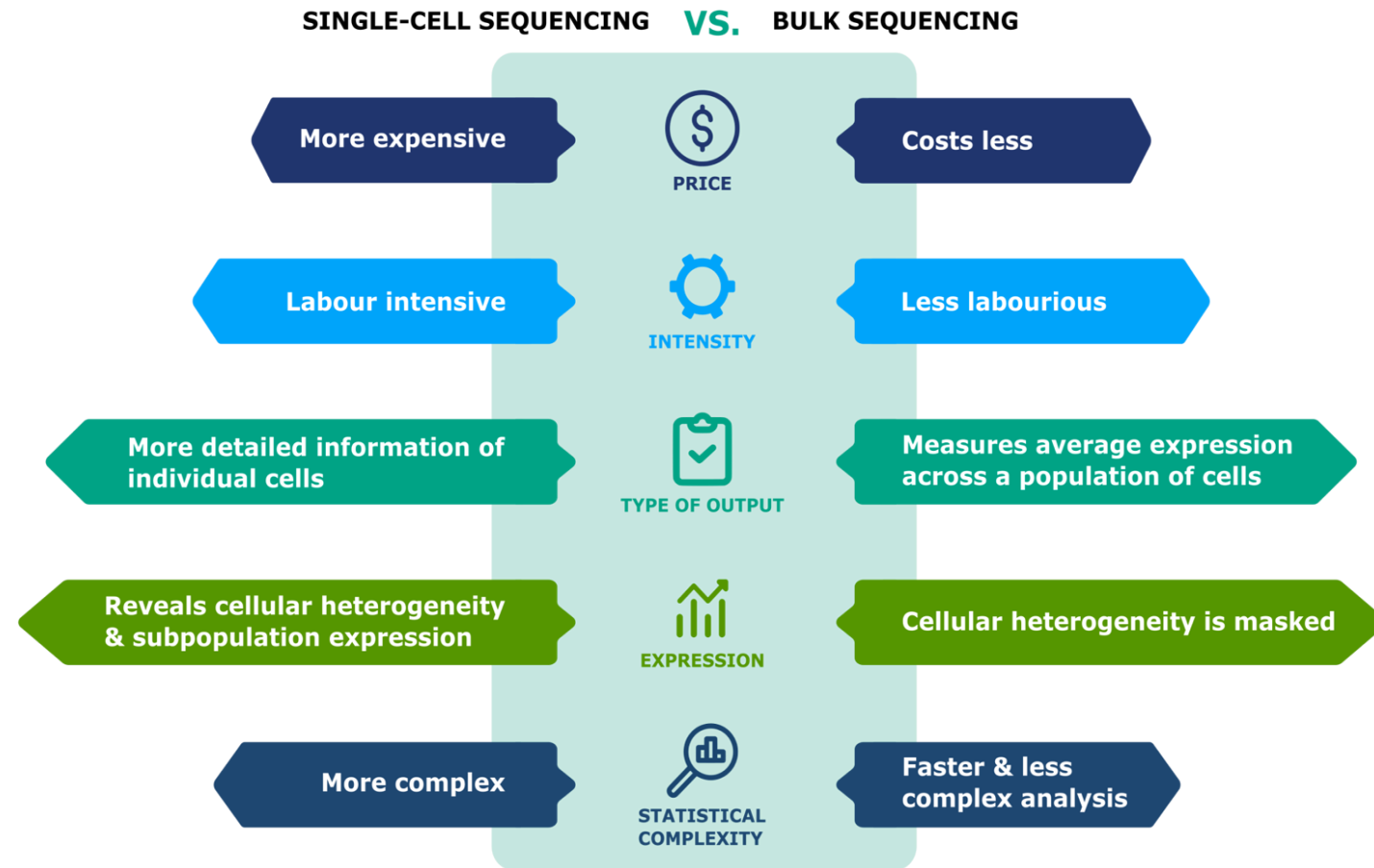

Challenges with scRNA-seq data

- Large volume of data
- Low depth of sequencing per cell (zero-inflation)
 - Often detecting only 10-50% of the transcriptome per cell

Increased complexity and richer datasets, means more room for misinterpretations and deriving wrong conclusions!

Recommendations

Do not perform single-cell RNA-seq unless it is necessary for the experimental question of interest.

→ Could you answer the question using bulk sequencing, which is simpler and less costly? Perhaps FACS sorting the samples could allow for bulk analysis?

Recommendations

Understand the details of the experimental question you wish to address.

→ Library preparation method and analysis workflow can vary based on the specific experiment and tissue

Recommendations

Avoid technical sources of variability, if possible:

Discuss experimental design with experts prior to the initiation of the experiment

Isolate RNA from all the samples at same time

Prepare libraries at same time or alternate sample groups to avoid batch confounding

Experimental design



Three fundamental pillars of experimental design



- Replication
- Randomization
- Blocking

Experimental design



Three fundamental pillars of experimental design

- Replication
- Randomization
- **Blocking**

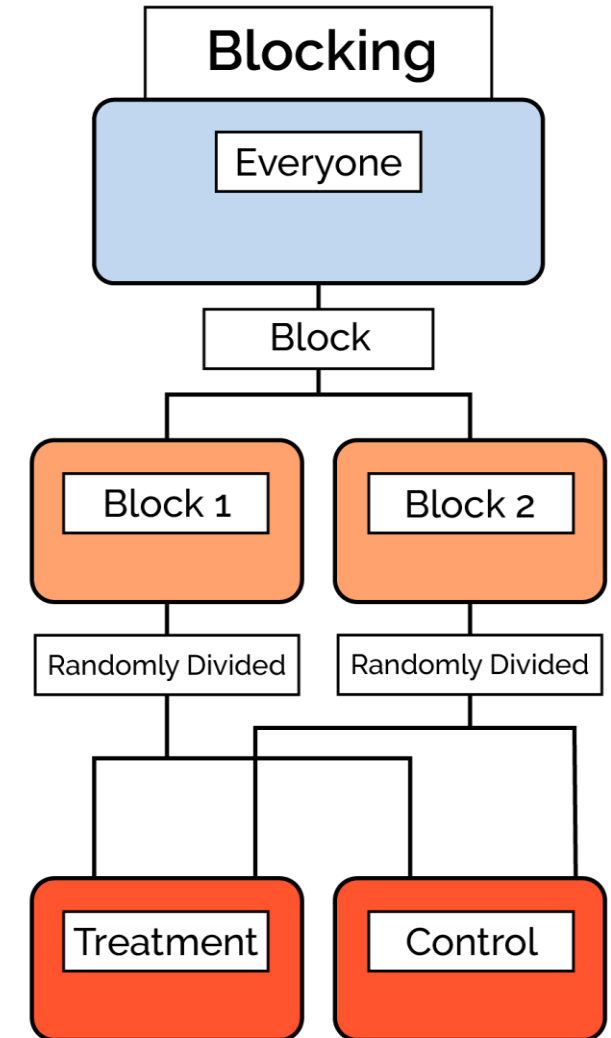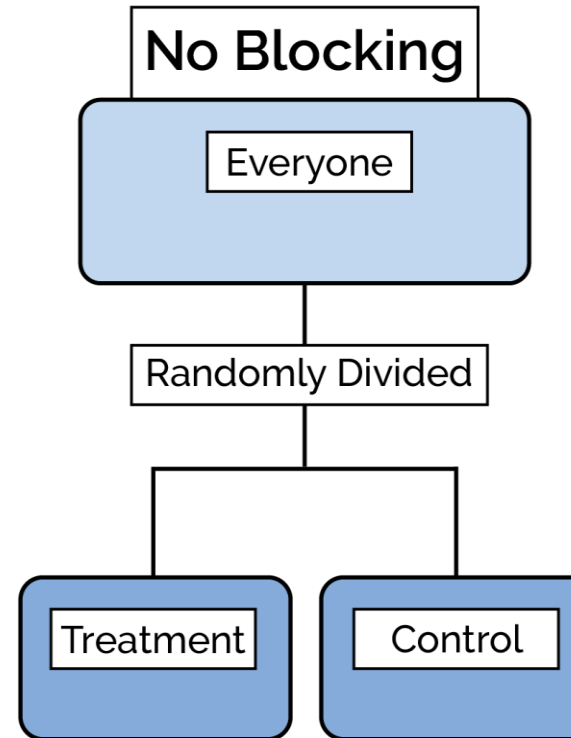

Experimental design



Replication, randomization and blocking

Be aware of confounding factors:

- » Person performing handling
- » Reagents
- » Sequencing lane/library

Further reading:

- » <https://doi.org/10.3389/fcell.2018.00108>
- » <https://doi.org/10.1093/bib/bby007>
- » <https://doi.org/10.1093/bfpg/elx035>

Experimental design



Replication, randomization and blocking

Be aware of confounding factors:

- » Person performing handling
- » Reagents
- » Sequencing lane/library

Record any factor for downstream correction:

Taking notes = the best QC

Further reading:

- » <https://doi.org/10.3389/fcell.2018.00108>
- » <https://doi.org/10.1093/bib/bby007>
- » <https://doi.org/10.1093/bfgp/elx035>

Experimental design



Replication, randomization and blocking

Be aware of confounding factors:

- » Person performing handling
- » Reagents
- » Sequencing lane/library

Record any factor for downstream correction:


Taking notes = the best QC

Different cells of one replicate $\neq$ replicates!

Further reading:

- » <https://doi.org/10.3389/fcell.2018.00108>
- » <https://doi.org/10.1093/bib/bby007>
- » <https://doi.org/10.1093/bfpg/elx035>

Sequencing output

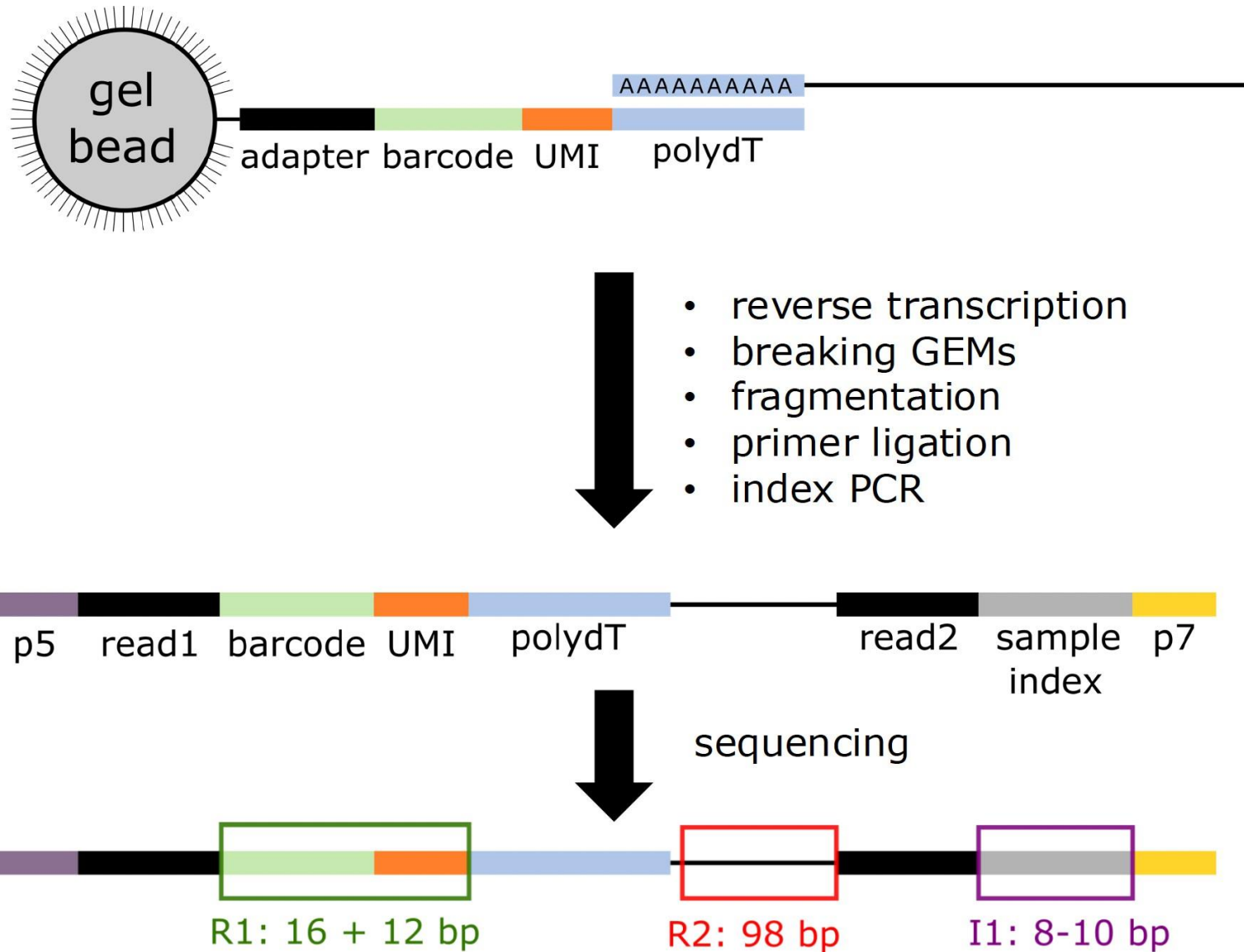

Sequencing output

sample#			read type	
ETV6-RUNX1_1_S1_L001_I1_001.fastq.gz				
ETV6-RUNX1_1_S1_L001_R1_001.fastq.gz				
ETV6-RUNX1_1_S1_L001_R2_001.fastq.gz				

sample ID lane

Index Read (I1) - contains cell barcodes + UMI

Read 1 (R1) - contains cell barcode + UMI (same as I1 in 10x v2/v3)

Read 2 (R2) - contains cDNA insert (gene/transcript sequence)

From FASTQ to counts

Tools for this part of the workflow include [Alevin](#), [UMI-tools](#), and [Cell Ranger](#) (10X data). While each tool will do things slightly differently the steps below are common to all:

1. Assign reads to cell
2. Alignment
3. Quantification: # UMI/gene • UMI: Unique Molecular Identifier
4. Cell calling

The 10X Software Suite

Pipeline for mapping, filtering, QC, and quantitation of libraries

Cell
Ranger

Barcode Extraction and filtering

- ⇒ Identifies cell level barcodes

Mapping to reference

- ⇒ Uses STAR aligner

Generate count table

- ⇒ UMIs per gene in each cell

Dimensionality Reduction

- ⇒ PCA and tSNE

Clustering

- ⇒ K-means and Graph Based

Cell Ranger

Cell Ranger Count (quantitates a single run)

```
$ cellranger count --id=COURSE \  
                    --transcriptome=/cellranger/references/GRCh38/ \  
                    --fastqs=/10X/ \  
                    --localcores=8 \  
                    --localmem=32
```

References:

- Human & mouse: download pre-built from 10x website
- Other organisms: custom reference with **cellranger mkref**

extensive documentation:

<https://support.10xgenomics.com/single-cell-gene-expression/software/pipelines/latest/what-is-cell-ranger>

Output files generated by Cell Ranger

`web_summary.html` - Web format QC report

`raw/filtered_features_bc_matrix`

- ❖ `barcodes.tsv.gz` - cell level barcodes seen in this sample
- ❖ `features.tsv.gz` - list of quantitated features (usually Ensembl genes)
- ❖ `matrix.mtx.gz` - (sparse) matrix of counts for cells and features

`possorted_genome_bam.bam` - BAM file of mapped reads

`molecule_info.h5` – Details of the cell barcodes

`cloupe.cloupe` - Analysis data for Loupe Cell browser

Cellranger report

ETV6-RUNX1_1

Alerts

The analysis detected 1 warning.

Alert	Value	Detail
 Fraction of RNA read bases with Q-score >= 30 is low	59.4%	Fraction of RNA read bases with Q-score >= 30 should be above 65%. A lower fraction might indicate poor sequencing quality. This is Read 1 for the Single Cell 3' v1 chemistry and Single Cell 5' paired end, Read 2 for the Single Cell 3' v2/v3 chemistry and Single Cell 5' R2-only)

Summary Analysis

3,091

Estimated Number of Cells

68,259

Mean Reads per Cell

1,717

Median Genes per Cell


Sequencing

Number of Reads	210,987,037
Number of Short Reads Skipped	0
Valid Barcodes	98.2%
Valid UMIs	100.0%
Sequencing Saturation	84.4%
Q30 Bases in Barcode	96.4%
Q30 Bases in RNA Read	59.4%
Q30 Bases in UMI	96.5%

Mapping

Reads Mapped to Genome	95.8%
Reads Mapped Confidently to Genome	92.9%
Reads Mapped Confidently to Intergenic Regions	5.2%
Reads Mapped Confidently to Intronic Regions	25.5%
Reads Mapped Confidently to Exonic Regions	62.2%
Reads Mapped Confidently to Transcriptome	58.2%
Reads Mapped Antisense to Gene	1.2%

Cells

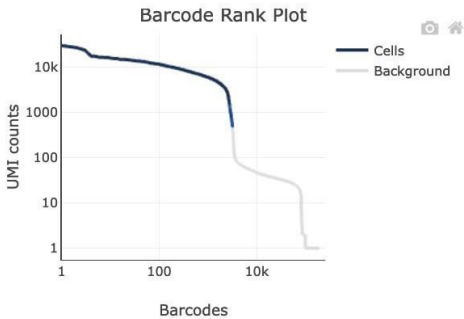


Estimated Number of Cells	3,091
Fraction Reads in Cells	91.1%
Mean Reads per Cell	68,259
Median Genes per Cell	1,717
Total Genes Detected	18,334
Median UMI Counts per Cell	4,825


Sample

Sample ID	ETV6-RUNX1_1
Sample Description	
Chemistry	Single Cell 3' v2
Include introns	False
Reference Path	...nger/refdata-cellranger-GRCh38-3.0.0
Transcriptome	GRCh38-3.0.0
Pipeline Version	cellranger-6.0.1



Cell calling

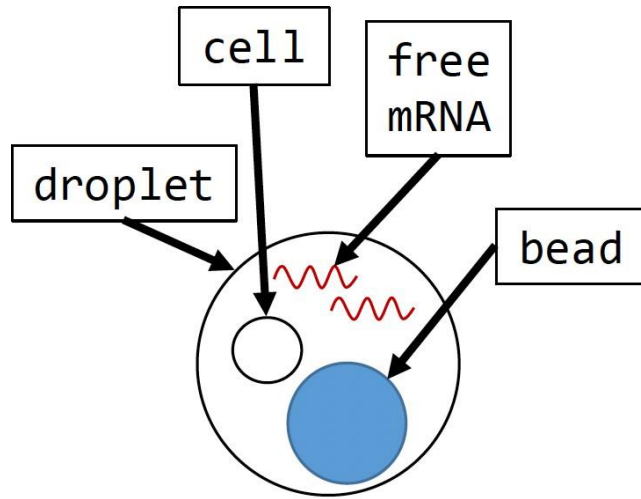

Cell
captured

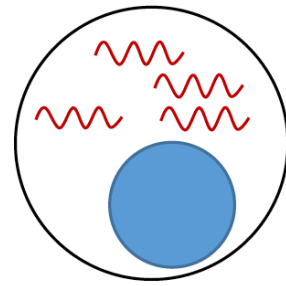

Only
background
captured

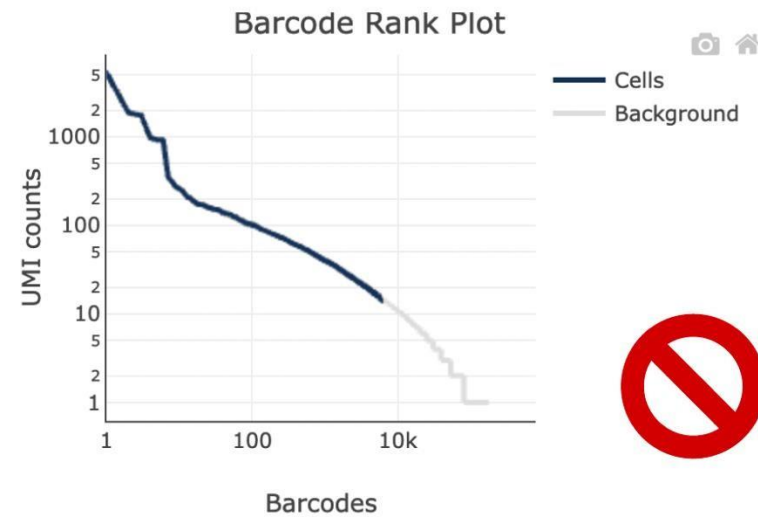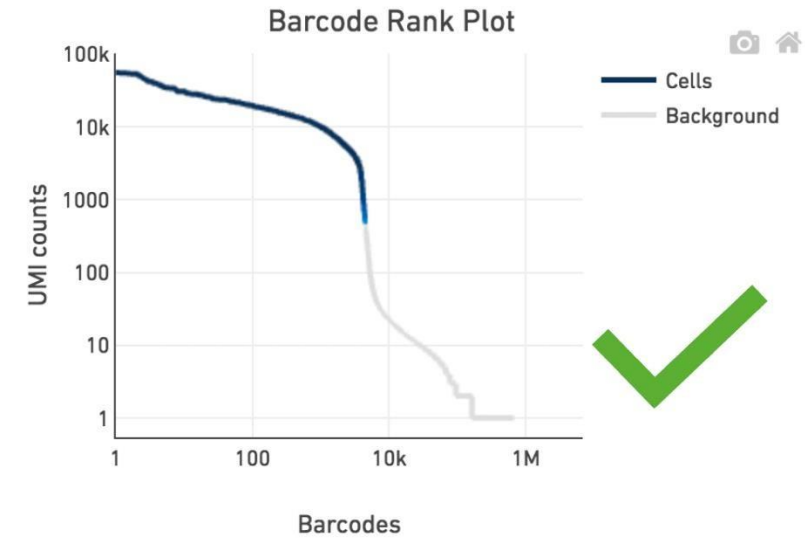

Background 'cells': low #UMI/barcode

Other important parameters

- Number of cells (<8% doublet rate):
 - <20k/channel (chromium X)
 - <10k/channel (chromium controller)


Other important parameters


- Number of cells (<8% doublet rate):
 - <20k/channel (chromium X)
 - <10k/channel (chromium controller)
- Reads per cell: 20k – 50k (depends on application)

Other important parameters

- Number of cells (<8% doublet rate):
 - <20k/channel (chromium X)
 - <10k/channel (chromium controller)
- Reads per cell: 20k – 50k (depends on application)
- Sequencing saturation

$$\text{saturation} = 1 - \frac{\# \text{ unique reads}}{\# \text{ reads}}$$

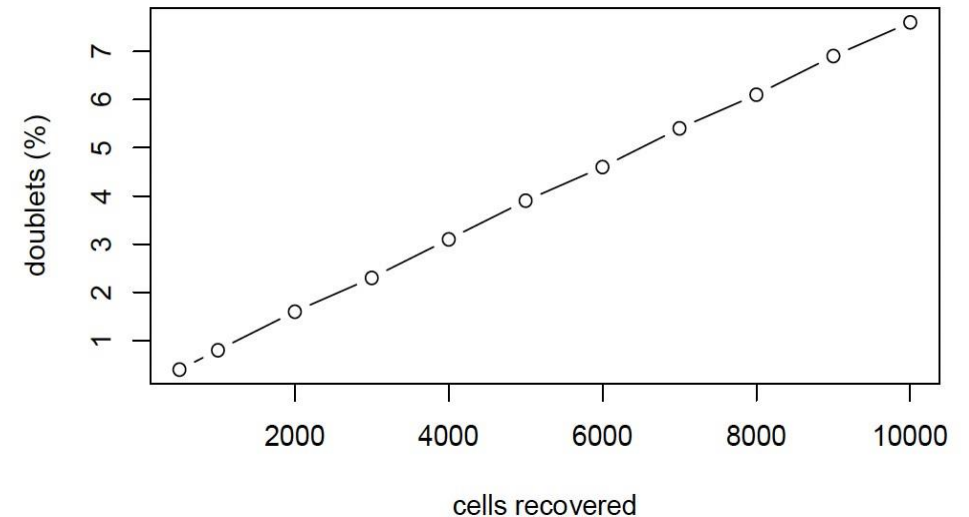

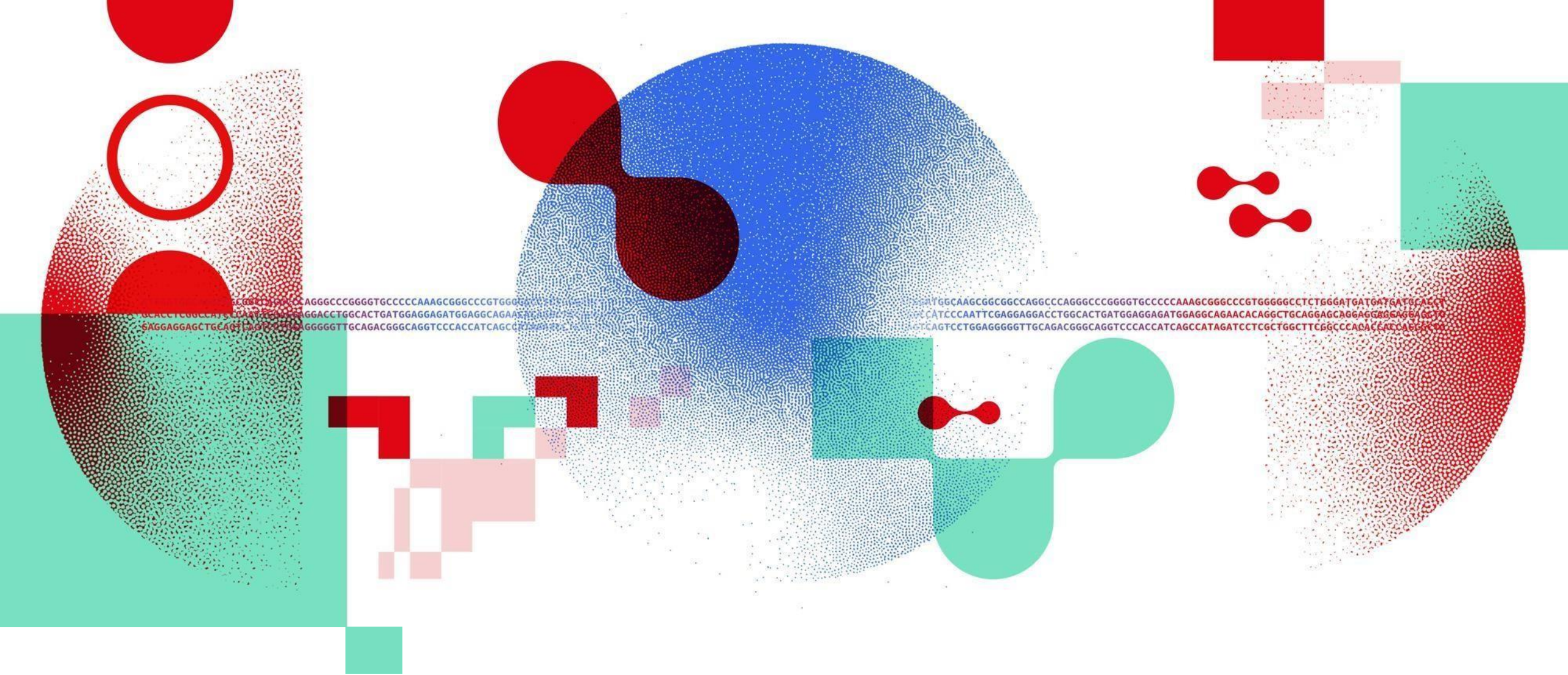

Thank you

DATA SCIENTISTS FOR LIFE

sib.swiss